

# Ambu® aBox™ 2

## DATASHEET

---

Ambu® aBox™ 2 is an intuitive all-in-one displaying and processing unit supporting your use of Ambu endoscopes.

You'll see a live endoscopic image within seconds when connecting an endoscope to aBox 2.

The intuitive user interface makes it easy to maintain focus and the adaptive light control and continuous image optimization allows you to see the details of interest.

aBox 2 has a built-in touchscreen providing you with a keyboard, access to image settings and showing the endoscopic image right there on the screen.

Videos and pictures can be exported via USB or directly to PACS via LAN or Wi-Fi.

With the option of importing a worklist the connection between patient and imagery is ensured when exporting the data.

aBox 2 can be used stand-alone or connected to available recorders and external screens.



## KEY BENEFITS

---

Full HD resolution LCD monitor

Adjustable touchscreen with a simple user interface

Continuous image optimization

Supports DICOM modality worklist

Compact and flexible design

Options to connect to external monitors and recorders

Built-in recorder

Easy to document imaging data

Export functions to a USB drive and PACS via DICOM

## TECHNICAL SPECIFICATIONS

### Display

Resolution	1920 x 1080 pixels
Orientation	Landscape
Display type	12.8" colour TFT LCD; capacitive touchscreen

### Connections

USB connections	Front: 1x USB 3.0; Back: 3x USB 3.0 and 2x USB 2.0 (all type A)
Video output	DVI-D (1920 x 1080p, 60 fps) and 3G-SDI (1920 x 1080p, 60 fps)
Wi-Fi	Supports IEEE Standards 802.11ac/a/b/g/n
LAN	Ethernet RJ45 connector, 10/100/1000 Mbps
Trigger output	2 x D-SUB9 and 2 x 3.5 mm Jack
Connection to potential equalization cable	DIN 42801 standard

### Memory

Storage Capacity	128 GB
------------------	--------

### Electrical power

Power requirement	100 - 240 VAC ; 2A
Protection against electrical shock	Class I

### Operating environment

Operating temperature	10 - 40 °C (50 - 104 °F)
Operating relative humidity	30 - 85 %
Operating atmospheric pressure	80 - 106 kPa
Operating altitude	≤ 2000 m

### Storage environment

Storage temperature	0 - 40 °C (32 - 104 °F)
Storage relative humidity	10 - 85 %
Storage atmospheric pressure	50 - 106 kPa

### IP Protection Classification

IP Protection Classification	IP31 - Protection against solid particles (3: < 2.5 mm) and liquid ingress (1: dripping water)
------------------------------	--

### Dimensions



Height - screen in 90° angle (a)	278 mm (10.9")
Height - screen folded flat (b)	121 mm (4.76")
Width (c)	330 mm (13.0")
Depth (d)	337 mm (13.3")
Weight	8 kg (17.6 lbs)

### Power cables

Power cables 6 interchangeable types	USA & Japan: Type B, Model NEMA, 5 AC grounded power plug Australian: Type I, AS3112, AC grounded power plug UK: Type G, BS1363, AC grounded power plug European: Type E/F, CEE 7, AC grounded power plug Danish: Type K, 2-5a, AC grounded power plug Swiss: Type J, AC grounded power plug
Length EU/UK/CH	1830mm (72") including connector
Length US/AUS/DK	1800mm (71") including connector

## ORDERING SPECIFICATIONS

Description	Reference number
aBox 2	505001000

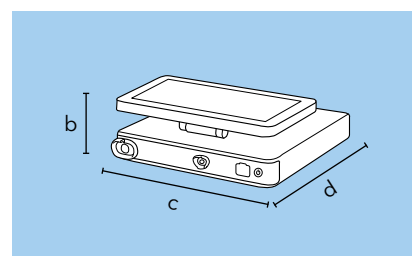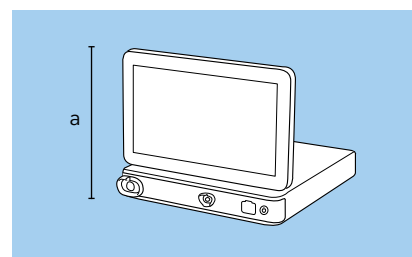

# Ambu

**Ambu A/S**  
Baltorpbakken 13  
2750 Ballerup  
Denmark  
T +45 72 25 20 00  
ambu.com



**UK RP Ambu Ltd**  
First Floor, Incubator 2  
Alconbury Weald Enterprise Campus  
Alconbury Weald  
Huntingdon PE28 4XA  
United Kingdom  
www.ambu.co.uk



Ambu is a registered trademark and aScope and aBox are trademarks of Ambu A/S.

US: Rx Only